Cystatin C Immunoparticles  
(ERM-DA471/IFCC Standardized)  
Code LX004

ENGLISH	
Intended use	For in vitro diagnostic use.  Cystatin C Immunoparticles (ERM-DA471/IFCC Standardized) are intended for the quantitative determination of cystatin C in human serum, heparinized plasma and EDTA plasma by turbidimetry. Cystatin C measurements are used as an aid in the diagnosis and treatment of renal diseases.
Summary and explanation	Cystatin C is a small, 13 kDa, non-glycosylated basic protein belonging to the cystatin super-family of cysteine protease inhibitors. Cystatin C is produced by virtually all nucleated cells, and is present in all investigated body fluids. The production rate is constant and is unaffected by inflammatory processes, gender, age and muscle mass (1). In the normal kidney, cystatin C is almost freely filtered through the glomerular membrane and then nearly completely reabsorbed and degraded by the proximal tubular cells. Therefore, the plasma concentration of cystatin C is almost exclusively determined by the glomerular filtration rate (GFR), making cystatin C an excellent indicator of GFR. Numerous studies and a meta-analysis incorporating 4492 subject samples have shown that serum cystatin C is superior to serum creatinine as a marker for GFR (2). In another meta-analysis incorporating 2007 subject samples, the diagnostic accuracy for impaired renal dysfunction was favored by serum cystatin C over serum creatinine (3).
Test principle	Human serum or plasma is mixed with the Cystatin C Immunoparticles. The resulting immune complexes are measured by turbidimetry. The signal generated is correlated with the concentration of cystatin C in the sample. By interpolation on a standard curve, the concentration of cystatin C in the sample is calculated.
Reagent	Purified immunoglobulin fraction of rabbit antiserum directed against cystatin C covalently coupled to uniform polystyrene particles. Provided as a suspension and preserved with 15 mmol/L sodium azide. The Cystatin C Immunoparticle suspension is ready for use.
Immunogen	Recombinant human cystatin C produced in <i>E. coli</i> was used as immunogen for raising the antibody coupled to the polystyrene particles.
Precautions	1. For in vitro diagnostic use. For professional users. 2. Interpretation of results must be made within the context of the patient's clinical history and other diagnostic tests by a certified professional. 3. The possibility of antigen excess, i.e. antigen concentrations exceeding the security range of the assay, should be considered. 4. This product contains sodium azide (NaN <sub>3</sub> ) a chemical highly toxic in pure form. At product concentrations, though not classified as hazardous, sodium azide may react with lead and copper plumbing to form highly explosive build-ups of metal azide. Upon disposal, flush with large volumes of water to prevent metal azide build-up in plumbing. 5. As with any product derived from biological sources, proper handling procedures should be used. 6. Wear appropriate Personal Protective Equipment to avoid contact with eyes and skin. 7. Unused solution should be disposed of according to local, State and Federal regulations.
Storage	Store at 2-8 °C. Do not use after expiration date stamped on vial. If reagents are stored under any conditions other than those specified, the user must verify the conditions. If unexpected results are obtained which cannot be explained by variations in laboratory procedures and a problem with the reagent is suspected, contact our Technical Services. For on-board stability data please refer to the Specific Performance Characteristics for the specific instrument.
Specimens	Sample material: Human serum, heparinized plasma and EDTA plasma. Samples can be stored at 2-8 °C, if analysed within 2 days. Otherwise, samples must be stored at -20 °C. Samples are stable for at least 6 months at -20 °C. It is recommended that frozen samples are thawed at 37 °C; samples must be mixed well before analysis.  As with any product derived from biological sources, proper handling procedures should be used. For interference please refer to the Specific Performance Characteristics for the specific instrument.
Procedure	<u>Materials provided</u> Dako Cystatin C Immunoparticles (ERM-DA471/IFCC Standardized) Code LX004  <u>Materials required but not provided</u> Dako Cystatin C Calibrator (ERM-DA471/IFCC Standardized) Code X7912 Dako Cystatin C Control Set (ERM-DA471/IFCC Standardized) Code X7913 Dako Reaction Buffer 9 Code S2361 Diluent See Application Note for the individual instrument  <u>Assay procedure</u> Detailed, validated Application Notes for the determination of cystatin C in human sample material, based on optimized test systems (4), are available from Dako for a range of automated instruments. For instruments not covered by a specific Application Note, a General Application Note, which provides general assay procedure guidelines is available. When using the General Application Note, each individual laboratory should validate their test system. Please consult the instruction manual of the instrument used regarding precautions, operation, maintenance, and instrument performance characteristics.

Establishment of the calibration curve

The assigned value of the Dako Cystatin C Calibrator (ERM-DA471/IFCC Standardized), Code X7912, is given on the Analytical Value Sheet included with the calibrator. The exact measuring range depends on the specific lot of the calibrator and the Application Note for the individual instrument.

Calibration curves are constructed by multi-point calibration. Depending on the instrument, a dilution series is either prepared automatically by the instrument or manually, and then applied on the instrument. For details on interval for recalibration, please refer to the Application Note the individual instrument.

Assay of samples

Samples are measured according to the Application Note for the individual instrument.

If the results obtained are outside the measuring range, the assay can be repeated using a post concentration or post dilution of the sample as stated in the Application Note.

Internal quality control

Dako Cystatin C (ERM-DA471/IFCC Standardized) Control 1 and Control 2, Code X7913, must be assayed with each establishment of a calibration curve and with each run of samples. The controls should be assayed and evaluated as samples. The assigned values and confidence intervals are given on the Analytical Value Sheet included with the controls. If quality control results fall outside the range stated on the Analytical Value Sheet or outside the range established by the individual laboratory, recalibration is recommended.

Results

The analytical results are automatically calculated by the instrument.

Limitations

For interference please refer to the Specific Performance Characteristics for the specific instrument.

Expected values

Reference range: 0.59-1.03 mg/L

The reference intervals were determined from a healthy population of 69 subjects ranging from 18 to 70 years in age. The reference intervals were calculated using non-parametrical statistics and represent the central 95% of the population. Healthy individuals ≤ 50 years of age have generally a glomerular filtration rate above 80 mL/min/1.73 m<sup>2</sup>, whereas the glomerular filtration rate of healthy individuals > 50 years decreases with age (5). The above reference intervals may be used as a guide for demonstrating normal or reduced GFR, but it is recommended to determine the reference intervals for the local population.

Specific performance characteristics

The Dako Cystatin C assay is based on the usage of Dako Cystatin C Immunoparticles together with Dako Cystatin C Calibrator (ERM-DA471/IFCC Standardized), Dako Cystatin C Control Set (ERM-DA471/IFCC Standardized), Dako Reaction Buffer 9, and the diluent recommended in the Application Note for the specific instrument. Please refer to the application notes for the performance characteristics of the Dako Cystatin C assay on specific instruments.

Specificity

The specificity of the cystatin C antibody, before coupling to the particles, was ascertained by crossed immunoelectrophoresis. Only the cystatin C arch appeared when the antibody was tested against concentrated human urine from patients with tubular proteinuria.

Method comparison

Results obtained with the Dako Cystatin C assay were compared with a commercially available turbidimetric assay.

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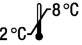
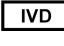
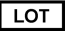
1. Grubb AO. Cystatin C – properties and use as diagnostic marker. Adv Clin Chem 2000;35:63-99.
2. Dharnidharka VR, Kwon C, Stevens G. Serum cystatin C is superior to serum creatinine as a marker of kidney function: A meta-analysis. Am J Kidney Dis 2002;40:221-6.
3. Roos JF, Doust J, Tett SE, Kirkpatrick CM: Diagnostic accuracy of cystatin C compared to serum creatinine for the estimation of renal dysfunction in adults and children-A meta-analysis. Clin Biochem 2007;40:383-91.
4. Blirup-Jensen S. Protein standardization III: Method optimization. Basic principles for quantitative determination of human serum proteins on automated instruments based on turbidimetry or nephelometry. Clin Chem Lab Med 2001;39:1098-1109.
5. Nilsson-Ehle P. Iohexol clearance for the determination of glomerular filtration rate: 15 years' experience in clinical practice. eJIFCC, Vol 13; No.2: <http://www.ifcc.org/ejifcc.htm>.

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